



Work instructions for AZ Calculation and Sampling

SiMPL EBR

V 1.0.0
February 2026



Table of Contents

1 PURPOSE	3
2 OVERVIEW	3
3 NAVIGATION	4
4 MASTER RECIPE EBR	5
1.1.1 Calculation	8
1.1.2 FIO/ Investigational Samples	10
1.1.3 DLIMS Samples	13
5 Revision History	15

1 Purpose

The purpose of this Work Instruction is to provide step-by-step guidance for executing the Calculation and Sampling process steps within the *Simple Manufacturing Process Layer Electronic Batch Records (SiMPL EBR)* application in G20. This instruction ensures accurate and consistent recording of preparation, verification, and usage data related to the specified process steps during manufacturing.

2 Overview


SiMPL EBR is a web-based application for batch record processing in medical product manufacturing. It integrates manufacturing processes and systems using SAP ECC's core functionality through the Fiori Launchpad.

This application combines pharmaceutical industry practices with the ECC framework, enabling flexible processing of EBRs. EBRs guide process operators through production steps and provide real-time process control visualization. Operators can maintain SAP EBRs, tracking process parameters, deviations, and digital signatures.

SiMPL EBR reduces the risks associated with outdated EBRs by running in any browser and utilizing standard process step and SAP integration capabilities.


3 Navigation

The *SiMPL EBR* application can be accessed through the SAP Fiori Launchpad in G20.

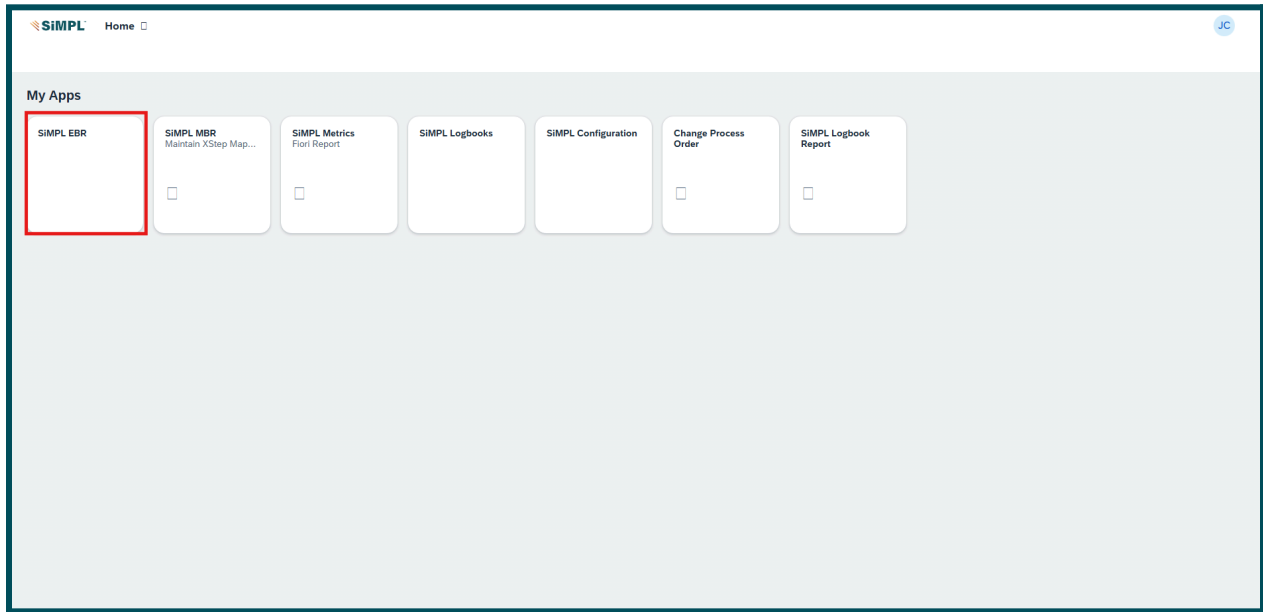

SiMPL EBR Application

The process orders are displayed only for the plant(s) to which the user has access. Other pages such as the Network Graph, Process Order Detail Page, and EBR Processing Page also work according to the user’s plant access. If a user does not have access to a specific plant, no process order details will be displayed on any page. Select the relevant Process Order on the Process Order List Page.

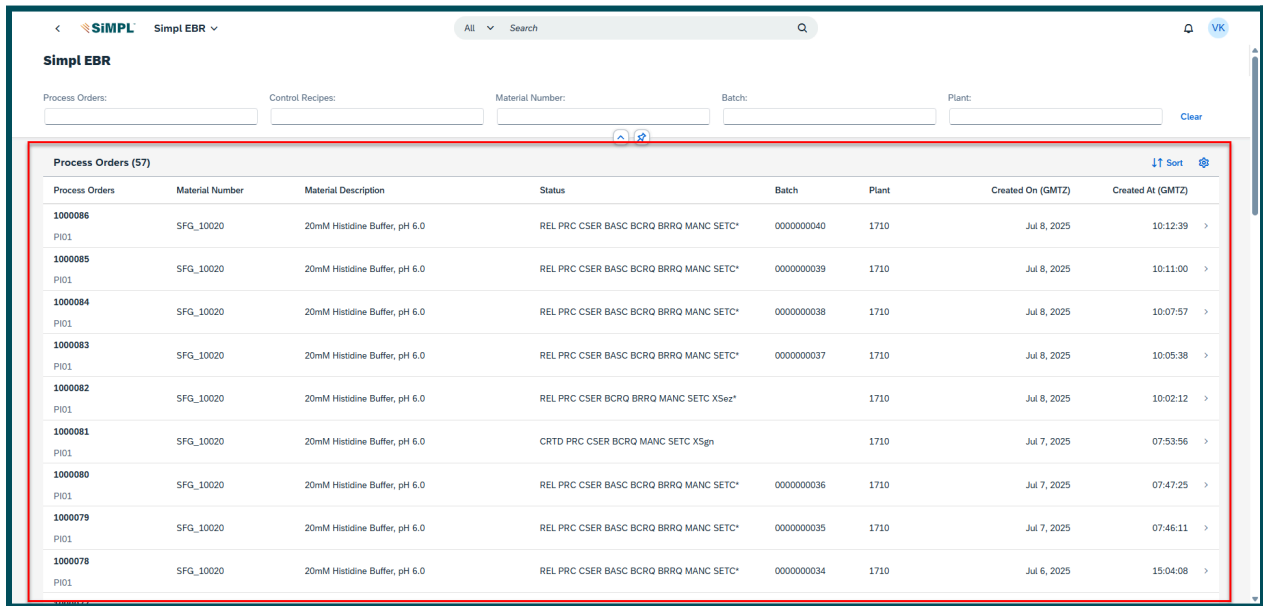

Process Orders	Material Number	Material Description	Status	Batch	Plant	Created On (GMTZ)	Created At (GMTZ)
1000086 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000040	1710	Jul 8, 2025	10:12:39 >
1000085 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000039	1710	Jul 8, 2025	10:11:00 >
1000084 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000038	1710	Jul 8, 2025	10:07:57 >
1000083 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000037	1710	Jul 8, 2025	10:05:38 >
1000082 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BCRQ BRRO MANC SETC X5ez*		1710	Jul 8, 2025	10:02:12 >
1000081 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	CRTD PRC CSER BCRQ MANC SETC X5gn		1710	Jul 7, 2025	07:53:56 >
1000080 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000036	1710	Jul 7, 2025	07:47:25 >
1000079 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000035	1710	Jul 7, 2025	07:46:11 >
1000078 PIQ1	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRRO MANC SETC*	0000000034	1710	Jul 6, 2025	15:04:08 >

Process Order List

4 Master Recipe EBR

To carry out production-related tasks, it is essential to access the correct process order in the FLP. This allows you to execute the EBR, which contains the step-by-step instructions required for processing. This section outlines the steps required to locate and open a specific process order within the Fiori Launchpad in order to initiate and execute the associated EBR efficiently.

1. Input the required process order into the search field and press Enter.
2. Once the relevant process order appears, select it to open and launch the associated EBR.

The screenshot shows the SIMPL Fiori Launchpad interface. At the top, there is a search bar with the text 'Simpl EBR'. Below the search bar, there are input fields for 'Process Orders', 'Control Recipes', 'Material Number', 'Batch', and 'Plant'. The 'Process Orders' field contains the value '1000083'. Below these fields, there is a table titled 'Process Orders (1)'. The table has columns for 'Process Orders', 'Material Number', 'Material Description', 'Status', 'Batch', 'Plant', 'Created On (GMTZ)', and 'Created At (GMTZ)'. The table contains one row with the following data: '1000083', 'SFG_10020', '20mM Histidine Buffer, pH 6.0', 'REL PRC CSER BASC BCRQ BRQ MANC SETC*', '0000000037', '1710', 'Jul 8, 2025', and '10:05:38'. A red box labeled '1' highlights the search input field, and a red box labeled '2' highlights the search results table.

Process Orders	Material Number	Material Description	Status	Batch	Plant	Created On (GMTZ)	Created At (GMTZ)
1000083 PI01	SFG_10020	20mM Histidine Buffer, pH 6.0	REL PRC CSER BASC BCRQ BRQ MANC SETC*	0000000037	1710	Jul 8, 2025	10:05:38

Steps 1-2: Process Order Details

3. Once the process order is loaded, the order details will be visible in the header section.
4. The radio buttons indicate whether the EBR is in 'Display' or 'Change' mode. If the operator intends to only view the EBR without performing any actions, it is recommended to keep it in 'Display' mode, as shown in the screenshot below. However, if the EBR will be executed, the operator must select the 'Change' mode option.
5. Click on the control recipe row to proceed to the associated EBR. This action opens the EBR interface, where users can either execute or review process instructions, depending on the selected mode (Display or Change).

Process Order
1000083

Order Information Control Recipes

Material: SFG_10020
Material Description: 20mM Histidine Buffer, pH 6.0
Plant: 1710
Version: 008

Recipe: 1
Recipe Group: 50000028
Document Type: P01
Batch: 0000000037

Total Quantity: 16.000
Delivered: 0.000
Storage Location:
Scheduled Start (GMTZ): Jul 9, 2025

Scheduled End (GMTZ): Jul 9, 2025
Status: REL PRC CSER BASC BCRO BRRO MANC SETC*
Process Order Type: Production

Control Recipes

Control Recipes (1) ☒ Display ☐ Change

Control Recipes	Description	Address	Status	Recipe Type	Control Recipe Destination	Destination Type	Created (GMTZ)	Changed On (GMTZ)
10000000000000000105	20mM Histidine Buffer, pH 6.0	MIXING	In process	No Test	MX	4	Jul 9, 2025 03:05:55	Jul 14, 2025

Steps 3-5: Display or Change mode

6. The EBR displays the process order with all phases.

SIMPL EBR Change Manufacturing

Process Order: 0000013612417
Material Number: 8008056
SIMPL Test MR

Control Recipe: 10000000000000364654
Batch: C8008


Phase 0020 Phase 0030 Phase 0040 Overview

- > 0020: Process Order Information
- > 0030: Weighout and Mixing
- > 0040: Room & Equipment Cleaning
- > Overview

Step 6: EBR


7. To expand the phases, click “Expand All” in the top-right corner.

8. To switch between “Display” and “Change” modes, click the **Change->Display** or **Display->Change** button.

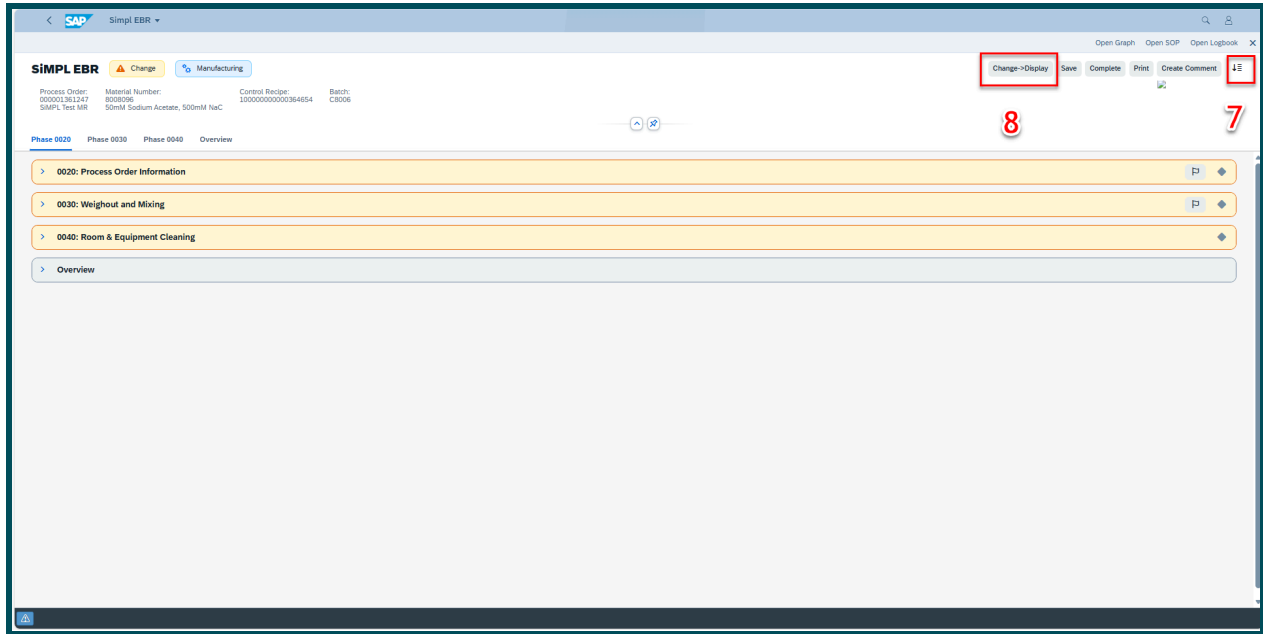

Steps 7-8: EBR

9. To collapse the phases view, click “Collapse All” in the top-right corner.

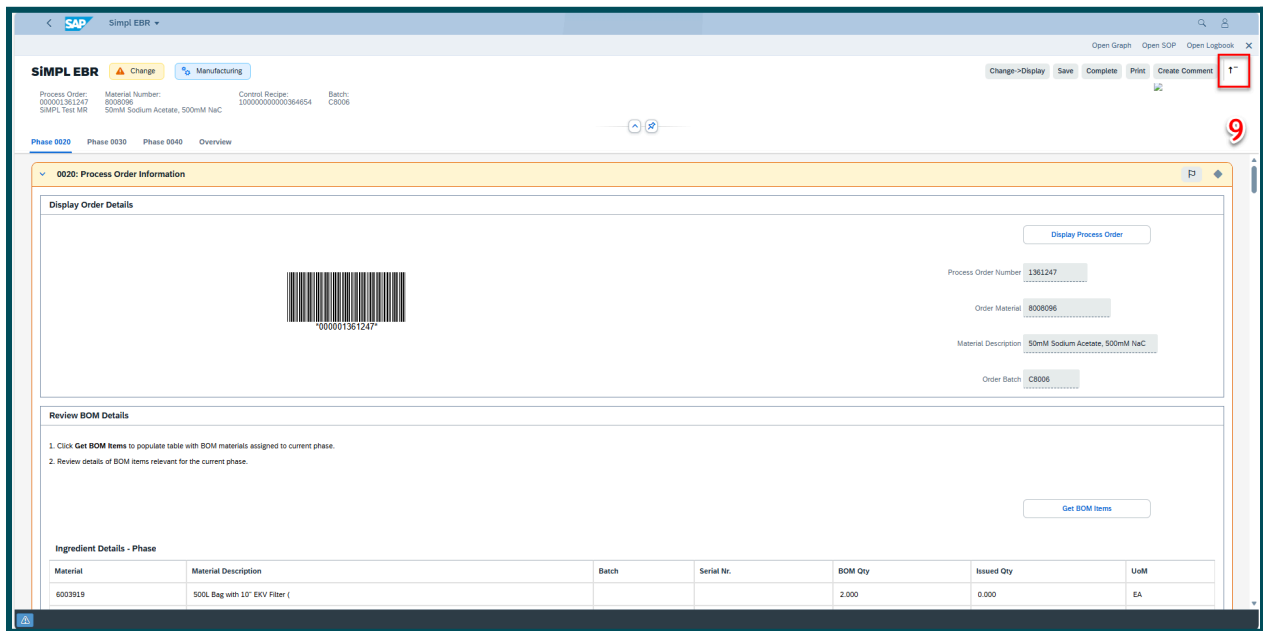

Step 9: Collapse All

1.1.1 Calculation

1. Enter the Measured Amount in the **Measured Amount** field and press enter.

The screenshot shows the SAP SIMPL EBR interface for the 'Calculation Polysorbate 20' form. The 'Measured Amount' field is highlighted with a red box and a red '1' above it, indicating where to enter the value. The 'Polysorbate 20 NF/EP/JPE Dens.' field contains the value '1.1100'. The 'Result' field contains the value '0'. The 'Performed By' and 'Witnessed By' fields are empty.

Measured Amount	Polysorbate 20 NF/EP/JPE Dens.	Result
0.000	1.1100	0

2. The results are calculated and populated in the **Result** field.

The screenshot shows the SAP SIMPL EBR interface for the 'Calculation Polysorbate 20' form. The 'Measured Amount' field now contains the value '20'. The 'Polysorbate 20 NF/EP/JPE Dens.' field contains the value '1.1100'. The 'Result' field is highlighted with a red box and contains the calculated value '18.018'. The 'Performed By' and 'Witnessed By' fields are empty.

Measured Amount	Polysorbate 20 NF/EP/JPE Dens.	Result
20	1.1100	18.018

3. After the user performing the tasks are satisfied with the entries, the user can sign the **Performed By** field.

The screenshot shows the SIMPL EBR interface for Phase 0030. The 'Calculation Polysorbate 20' section contains a table with the following data:

Measured Amount		Polysorbate 20 NF/EP/LPE Dens.	Result
20	/	1.1100	18.018

Below the table, the 'Performed By' field is highlighted with a red box and a red circle with the number 3. The 'Witnessed By' field is also visible below it.

4. The supervisor responsible for confirming the entry and its verification must sign the **Witnessed By** field, ensuring that this action is performed by a different user than the one who signed **Performed By**.

The screenshot shows the SIMPL EBR interface for Phase 0030. The 'Calculation Polysorbate 20' section contains a table with the following data:

Measured Amount		Polysorbate 20 NF/EP/LPE Dens.	Result
20	/	1.1100	18.018

Below the table, the 'Witnessed By' field is highlighted with a red box and a red circle with the number 4. The 'Performed By' field is also visible above it.

1.1.2 FIO/ Investigational Samples

1. Click the **Get Equipment** button to select an assigned BSC.

The screenshot shows the SIMPL EBR interface for Phase 0060. The 'Sampling' section contains a text box with instructions: 'These samples are for Investigational purposes only. The samples will only be tested for sterility if an investigation is initiated for the cell bank. Notify CCFS-TTCC for storage and testing.' Below this, the 'Get Equipment' button is highlighted with a red box, and a red '1' is placed next to it. Below the button are two input fields: 'Equip ID' and 'Equip Description'. At the bottom, there is a table with columns: Test, Sample Volume, Tube Type, Sample Pulled, Storage Temperature, Physical Storage Location, and Performed By. The table has one row for 'Sterility' with values: 1 x 20 ml, 50 mL Centrifuge tube, a dropdown menu, 2°C - 8°C, and empty fields for Physical Storage Location and Performed By.

Test	Sample Volume	Tube Type	Sample Pulled	Storage Temperature	Physical Storage Location	Performed By
Sterility	1 x 20 ml	50 mL Centrifuge tube		2°C - 8°C		

2. Once selected, the equipment selected should be retrieved and pulled into the **Equip ID** and **Equip Description** fields.

The screenshot shows the SIMPL EBR interface for Phase 0060. The 'Sampling' section contains the same text box as the previous screenshot. Below it, the 'Equip ID' field is populated with '30021754' and the 'Equip Description' field is populated with 'BIOSAFETY CABINET M205844 B...'. These two fields are highlighted with a red box, and a red '2' is placed next to it. The 'Get Equipment' button is still visible above the fields. The table at the bottom remains the same as in the previous screenshot.

Test	Sample Volume	Tube Type	Sample Pulled	Storage Temperature	Physical Storage Location	Performed By
Sterility	1 x 20 ml	50 mL Centrifuge tube		2°C - 8°C		

3. Select an option from the **Sampled Pulled** dropdown menu.

The screenshot shows the SIMPL EBR interface with the Phase 0060 tab selected. The Sampling table is visible, and the 'Sample Pulled' dropdown menu is highlighted with a red box and a red number 3. The table contains the following data:

Test	Sample Volume	Tube Type	Sample Pulled	Storage Temperature	Physical Storage Location	Performed By
Sterility	1 x 20 ml	50 mL Centrifuge tube	Yes	2°C - 8°C	Prep Room	
Retain Sterility	1 x 20ml	50 mL Centrifuge tube	Yes	2°C - 8°C	Prep Room	

4. Enter the **Storage Location** field.

The screenshot shows the SIMPL EBR interface with the Phase 0060 tab selected. The Sampling table is visible, and the 'Physical Storage Location' field is highlighted with a red box and a red number 4. The table contains the following data:

Test	Sample Volume	Tube Type	Sample Pulled	Storage Temperature	Physical Storage Location	Performed By
Sterility	1 x 20 ml	50 mL Centrifuge tube	Yes	2°C - 8°C	Prep Room	
Retain Sterility	1 x 20ml	50 mL Centrifuge tube	Yes	2°C - 8°C	Prep Room	

- Once the recorded data is reviewed and completed, the performer must sign in the relevant **Performed By** field.

The screenshot shows the SIMPL EBR SAP interface. The top navigation bar includes the SAP logo, 'Simpl EBR', and search/user icons. Below the navigation bar, there are tabs for 'Phase 0020', 'Phase 0030', 'Phase 0040', 'Phase 0050', 'Phase 0060' (selected), 'Phase 0070', and 'Overview'. The main content area displays a 'Get Equipment' button, 'Equip ID' (30021754), and 'Equip Description' (N/A). Below this is a 'Sampling' table with two rows: 'Sterility' and 'Retain Sterility'. The 'Performed By' field in the 'Sterility' row is highlighted with a red box and a red number 5.

Test	Sample Volume	Tube Type	Sample Pulled	Storage Temperature	Physical Storage Location	Performed By
Sterility	1 x 20 ml	50 mL Centrifuge tube	Yes	2°C - 8°C	Prep Room	
Retain Sterility	1 x 20ml	50 mL Centrifuge tube	Yes	2°C - 8°C	Prep Room	

1.1.3 DLIMS Samples

1. Press the **Get # of the Vessels** button.

The screenshot shows the SIMPL EBR interface. At the top, there's a header with 'SAP' and 'Simpl EBR'. Below it, there are tabs for 'Phase 0020', 'Phase 0030' (selected), 'Phase 0040', and 'Overview'. The main content area is titled 'Sampling & Test Results'. It contains instructions: 'Submit a 10 mL sample of this solution from each vessel in a PET container(s) for ENDOTOXIN testing. Use DLIMS form type "G".', 'PET= Polyethylene Terephthalate', and 'G = Green Form for development sample testing'. A red box with a red number '1' highlights the 'Get # of Vessels' button. Below this is a 'Sample Table' with columns: 'Vessel #', 'DLIMS Project Number', 'DLIMS Sample Number', and 'Sampled By'. The table has two rows: Row 1 with Vessel # 1, DLIMS Project Number 5000000268, and DLIMS Sample Number 5000000368; Row 2 with Vessel # 2, DLIMS Project Number 5000000267, and DLIMS Sample Number 5000000367. The 'Sampled By' column is empty. At the bottom right, there is a 'Witnessed By' field.

Vessel #	DLIMS Project Number	DLIMS Sample Number	Sampled By
1	5000000268	5000000368	
2	5000000267	5000000367	

2. Enter the **DLIMS Project Number** and **DLIMS Sample Number** fields.

The screenshot shows the same SIMPL EBR interface as before. The 'Sample Table' is now highlighted with a red box and a red number '2'. The table contains the same data as in the previous screenshot. The 'Get # of Vessels' button is still visible but not highlighted.

Vessel #	DLIMS Project Number	DLIMS Sample Number	Sampled By
1	5000000268	5000000368	
2	5000000267	5000000367	

3. Once the recorded data is reviewed and completed, the performer must sign in the relevant **Sampled By** field.

SAP Simpl EBR

Phase 0020 Phase 0030 Phase 0040 Overview

SIMPL EBR Change Manufacturing

Change->Display Save Complete Print Create Comment

Open Graph Open SOP Open Logbook

Sampling & Test Results

Submit a 10 mL sample of this solution from each vessel in a PET container(s) for ENDOTOXIN testing. Use DLIMS form type "G".
 PET= Polyethylene Terephthalate
 G = Green Form for development sample testing

Get # of Vessels

Sample Table

Vessel #	DLIMS Project Number	DLIMS Sample Number	Sampled By
1	5000000268	5000000368	
2	5000000267	5000000367	

Witnessed By

- The supervisor responsible for confirming the entry and its verification must sign the **Witnessed By** field, ensuring that this action is performed by a different operator than the one who signed **Performed By**. Once the **Witnessed By** field is signed, the process step is considered complete, and no further modifications can be made.

SAP Simpl EBR

Phase 0020 Phase 0030 Phase 0040 Overview

SIMPL EBR Change Manufacturing

Change->Display Save Complete Print Create Comment

Open Graph Open SOP Open Logbook

Sampling & Test Results

Submit a 10 mL sample of this solution from each vessel in a PET container(s) for ENDOTOXIN testing. Use DLIMS form type "G".
 PET= Polyethylene Terephthalate
 G = Green Form for development sample testing

Get # of Vessels

Sample Table

Vessel #	DLIMS Project Number	DLIMS Sample Number	Sampled By
1	5000000268	5000000368	
2	5000000267	5000000367	

Witnessed By

Craig Jason

Sampling

Test	Sample Volume	Tube Type	Sample Pulled	Storage Temperature	Storage Location	Performed By

5 REVISION HISTORY

Version No.	Description	Author	Date
1	Initial Document	Ruby Lutwama	14/05/2026